Body Fat Analyzer

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BFA-90

OWNER'S MANUAL

KULLANIM KILAVUZU

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Thank you very much for selecting the CAS Body Fat Analyzer BFA-90.

Please do read the user manual carefully and thoroughly so as to ensure the safe usage of this product, and keep the manual well for further reference in case you have problems.

TABLE OF CONTENT

Frequently asked questions.....	12
---------------------------------	----

Schedule for reference

- Body fat percentage	12
- Muscle percentage	12
- Bone mass	12
- Body Water (Hydration) percentage	13
- BMI (body mass index)	13

Specification	13
---------------------	----

Key function.....	14
-------------------	----

How to use	14
------------------	----

Manual	16
--------------	----

Attention	17
-----------------	----

Maintenance	18
-------------------	----

Warranty	18
----------------	----

Frequently asked questions

Is there anything wrong with "----" displayed on the screen?

- The scale is auto zero setting, please wait for a second, "0" will be displayed.
- No parameter was found in scale, please input your personal parameter.

The scale is only weight displayed, no any body fat percentage is found?

When measuring the body fat, shaking will disturb the body fat percentage measuring.
Please keep your body straight when weighing, and be sure that your foot bare and keep in good contact with metal sheet on scale.

The body fat percentage is much higher than normal?

Please double-check your personal parameter covering gender, height and age, make sure your bare foot is clear and dry, have a good contact with metal sheet, and big loss of water, quite exciting condition or callus on foot will affect the reading of body fat scale.

Schedule for reference

Body fat percentage (%)

Age	Female				Male			
	Underfat (-)	Healthy (0)	Overfat (+)	Obese (++)	Underfat (-)	Healthy (0)	Overfat (+)	Obese (++)
10~20	≤19.5	19.6~33.5	33.6~38.5	≥38.6	≤9.5	9.6~20.5	20.6~25.5	≥25.6
21~30	≤20.0	20.1~34.0	34.1~39.0	≥39.1	≤10.0	10.1~21.0	21.1~26.0	≥26.1
31~40	≤20.5	20.6~34.5	34.6~39.5	≥39.6	≤10.5	10.6~21.5	21.6~26.5	≥26.6
41~50	≤21.0	21.1~35.0	35.1~40.0	≥40.1	≤11.0	11.1~22.0	22.1~27.0	≥27.1
51~60	≤21.5	21.6~35.5	35.6~40.5	≥40.6	≤12.0	12.1~23.0	23.1~28.0	≥28.1
61~70	≤22.0	22.1~36.0	36.1~41.0	≥41.1	≤13.0	13.1~24.0	24.1~29.0	≥29.1
71~80	≤22.5	22.6~36.5	36.6~41.5	≥41.6	≤14.0	14.1~25.0	25.1~30.0	≥30.1

Muscle percentage (%)

Muscle	Female	Male
	>30	>35

Bone mass (kg)

	Female			Male		
	<45	45~60	>60	<60	60~75	>75
bone mass	1.8	2.2	2.5	2.5	2.9	>3.2

Body Water (Hydration) percentage (%)

Age	Female				Male			
	Low (-)	Normal (0)	Overhigh (+)	High (++)	Low (-)	Normal (0)	Overhigh (+)	High (++)
10~20	≤56.0	56.1~62.0	62.1~68.0	≥68.1	≤58.0	58.1~65.0	65.1~72.0	≥72.1
21~30	≤54.0	54.1~60.0	60.1~66.0	≥66.1	≤56.0	56.1~63.0	63.1~70.0	≥70.1
31~40	≤51.0	51.1~57.0	57.1~63.0	≥63.1	≤54.0	54.1~61.0	61.1~68.0	≥68.1
41~50	≤48.0	48.1~54.0	54.1~60.0	≥60.1	≤51.0	51.1~58.0	58.1~65.0	≥65.1
51~60	≤45.0	45.1~51.0	51.1~57.0	≥57.1	≤48.0	48.1~55.0	55.1~62.0	≥62.1
61~70	≤42.0	42.1~48.0	48.1~54.0	≥54.1	≤45.0	45.1~52.0	52.1~59.0	≥59.1
71~80	≤38.0	38.1~44.0	44.1~50.0	≥50.1	≤42.0	42.1~49.0	49.1~56.0	≥56.1


BMI (body mass index) = weight (kg) / (height)² (m²)%

BODY	Underfat	Standard	Overfat	Obese
BMI	< 16.5	16.5~25	25~30	> 30

Specification

Product Name	Body Fat Analyzer (BFA-90)
Equipped with a high precision strain gauge sensor system	
Max. Weight	180kg
Min. Effective	8kg
Units	KG
BMI/Body fat/Muscle/Water division	0.1%
Body bone division	0.1kg
Age range for body fat	10~80 years
Height range	100~220cm/3ft-3in~7ft-3in
Body bone range	0~15kg
Muscle range	13.0~100.0%
Fat range	5~50%
Water range	35~75%
Height unit	cm/ft
Memory	4 persons
Over-load / Low power indication	
Auto zero resetting / Auto power off	
Power	Automatic power generation control
Battery-free technology	


Key function

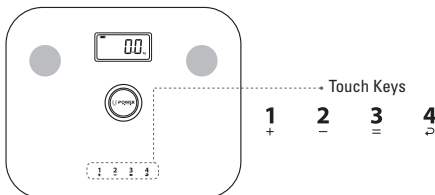
- 1** + to increase the data or swift the gender.
- 2** - to reduce the data or swift the gender.
- 3** = to confirm the data input.
- 4** ↻ to go back to last step in case any mistake from above operation.mistake from above operation.
Press one of the buttons for a long time if you want to adjust the parameter quickly.

How to use

Before using, put the scale on the hard and flat surface. There are two modes to start.

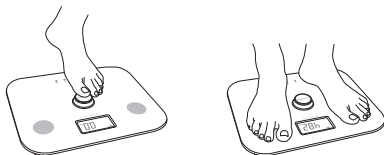Mode 1: For weighing only,no ingredient analysis

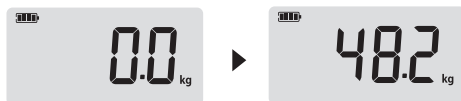- At least two quick stamping on the Magic-Start button will
- Generate power for weighing,please wait until "0.0" appears.

! Do not step on the scale until "0.0" appears.

Mode 2: For weighing and ingredient analysis

for weighing and ingredient analysis.

Step 1: At least two quickly stamping on the power button by your feet to create energy Family user number selection, press any of the key **1** **2** **3** **4** for a long time to setting and input your related personal data into the memory before using for the first time.


Step 2: Gender Selection

When the "MALE" is twinkling

Press **1** or **2** to select male or female

Press **3** to confirm

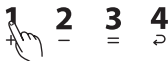
Step 3: Units setting :

When the unit is twinkling

Press **1** or **2** to change unit

Press **3** to confirm


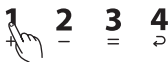
Step 4: Height setting :

When the height number is twinkling

Press **1** or **2** to change height

Press **3** to confirm

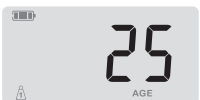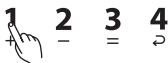
Step 5: Age setting :

When the age number is twinkling

Press **1** or **2** to change age

Press **3** to confirm

Press one of the buttons for a long time if you want to adjust the parameter quickly.

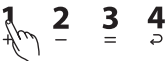Step 6 : Weighing and measuring:

When "0..." is displayed from left side to right side. It is analyzing.
Please do not leave the scale until your weight is locked on the screen.

NOTE

- 1)When "0.0" appears, stand with clean and dry bare feet on the scale, Your feet must stand evenly on the electrodes and do not stand outside of the electrodes.
- 2)Begin from the second using, you can press **1 2 3 4** for few seconds to begin measuring directly
- 3)This button "UNIT" can not change the units, please change the units on the Step 2 of mode 2.

Manual

Error Indication

Over-load Indication

When the scale is over-loaded[more than 150kg], the screen will appear "Err1".

It means the wrong operation caused the fat percentage "Err2".

Low battery indication

When the screen indication "  ", it will turn off automatically after few seconds.
Simply reload it by stamping on the U-power button!

Attention

Attention for safety

- The scale could not be used for the person have a pacemaker inside or with any other medical transplanting. The impedance passed through our body will have disturbance of the normal working of medical equipment.
- This scale is for family use only, could not do the medical job with any other medical equipment, we do not advise to use it as a medical equipment.
- Do not use it on the wet and slippery floor, do not stand on at the edge of scale to avoid any falling down.
- Advise disable person, pregnant woman, old person to use the scale with help.

Attention for operation

The body fat reading will be varied with the amount of body water and following factor will affected body fat percentage reading, diet style, drinking too much alcohol, heavy exercise, certain prescription drugs or diuretics, illness, or woman's menstrual cycle.

- Please take off your socks and keep your bare foot clean and dry.
- Advice to use the scale without any wearing, your wearing will influence the accuracy of the reading.
- Over diet or water lossing will have negative effect on the reading.
- Heavy exercise will have negative effect on the reading.
- Unstable or incorrect reading is aroused by not good connect with the metal sheet, please clean the dirt on the scale and step on scale as to instruction.
- Keep stable when weighing, any shaking will lead to the unaccuracy of scale weighing.
The scale could not be applied to following groups:
 1. People with age below 10 years old or above 80 years old.
 2. Athletic or muscular body type(body builder)
 3. Pregnant woman, person with dropsy.
- For above group, the reading of scale is only for reference.

REMARKS

Due to continuous improvements, this product is subject to change without notice.
Product appearance, color, if changes, according to the real object please.

Maintenance

Please place the scale on a clean and level surface.

- Please do not use the scale in extreme high or low temperature environment.
- Please do not jump or stamp on the scale and avoid dropping it, as the scale is a highly precise electronic instrument.
- To prevent scratches on glass panel, it is advised not to wear any footwear while taking any measurements.
- When cleaning the scale, please use only normal cloth and no chemical cleaners should be added to avoid any damage to the scale.
- This scale is only designed for home use and should not be used commercially for any reason at all. The < Body Fat Ratio Chart > is for reference purpose only.
- Please remove battery from the scale if it is not used frequently.

Warranty

- CAS warrants its products free of defects in materials and workmanship in normal use for a period of ONE year from the date of retail purchase.
- This warranty does NOT cover damages caused by misuse or abuse, including but not limited to:
 - Failure caused by unauthorized repairs or modifications;
 - Damage caused by shock or drop during transportation;
 - Failure caused by improper operation inconsistent with the instructions stated in this user manual;
 - Malfunction or damage from failure to provide the recommended maintenance;
 - Damage caused by improper use of power supply.
- Should this device require maintenance (or replacement at our option) under warranty, please deliver the original package to CAS prepaid. Please return the store receipt (with the retail purchase date) and a note with reasons to return on it as well.

Body Fat Analyzer

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BFA-90

Customer Support Center

+ 90 216 540 81 20

CAS Elektronik San. ve Tic. A.Ş.
Şerifali Mah. | Mevdudi Sok.No:26 |
Y.Dudullu – Ümraniye / İstanbul | 34775 | Türkiye
Tel. + 90 216 540 81 20

CAS